


Root[®] with EMMA[™] Bluetooth[®]

Seamless Integration of Tetherless Mainstream Capnography




Tetherless Solution

Cable-free design and quick Bluetooth pairing facilitate easy movement between care areas, for all patient populations

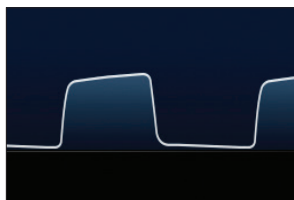

Increased Visibility

High-resolution, secondary display of large, crisp waveforms and numbers improves visibility of waveform for better interpretation

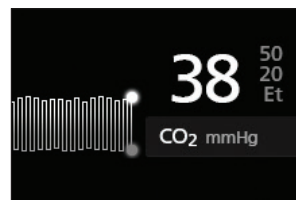

Real-time Trending

Clear trend display helps clinicians review progress over time and guide ventilation efforts



Automated EMR Documentation

Electronic charting of patient data in the EMR via Masimo Patient SafetyNet[™]* or Iris Gateway[™] simplifies workflow

Use of Capnography for Quick Patient Assessment


American Heart Association guidelines recommend continuous quantitative waveform capnography in addition to clinical assessment to help clinicians:



Confirm endotracheal tube placement¹




Assess the depth and effectiveness of chest compressions during CPR²

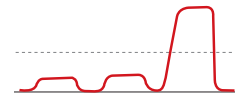

Recognise the return of spontaneous circulation during CPR²

Supplemental Capnography Display on Root

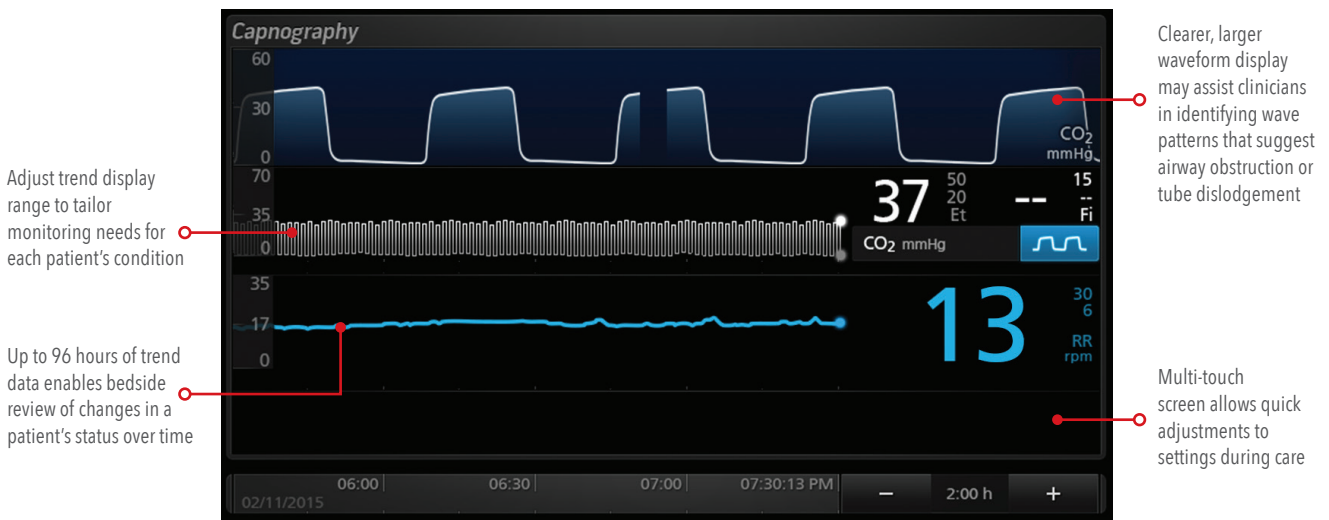

Quick assessment in a dedicated, adjustable window



Audible and visual alarms for quick identification of alarming parameters



One-touch Bluetooth menu access for pairing or re-pairing to multiple Root monitors

Pairing in Three Simple Steps

- 1 From the Root Bluetooth screen, select the EMMA Pair button
- 2 On EMMA, press power button twice until the screen displays the Bluetooth symbol
- 3 Select EMMA on the Devices Found list on Root and touch the Pair button

Ordering Information

ROOT	EMMA BLUETOOTH
Root9515	EMMA Bluetooth (mmHg)4271
Root with noninvasive blood pressure and temperature9695	EMMA Bluetooth (kPa).....4270
	EMMA Airway Adapter Adult/Paediatric17448
	EMMA Airway Adapter Infant.....17449

¹ Link MS, et al. *Circulation*. 2015; 132(suppl 2): S444-S464. ² Neumar RW et al. *Circulation*. 2010;122:S729-S767.
 * The use of the trademark Patient SafetyNet is under license from University Health System Consortium.

EMMA Bluetooth is not licensed for sale in Canada.

For professional use. See instructions for use for full prescribing information, including indications, contraindications, warnings, and precautions.

Masimo U.S.
 Tel: 1 877 4 Masimo
 info-america@masimo.com

Masimo International
 Tel: +41 32 720 1111
 info-international@masimo.com

